


medartis®

PRECISION IN FIXATION

PRODUCT INFORMATION

# Radial Head Plates 2.0

APTUS®  
Elbow


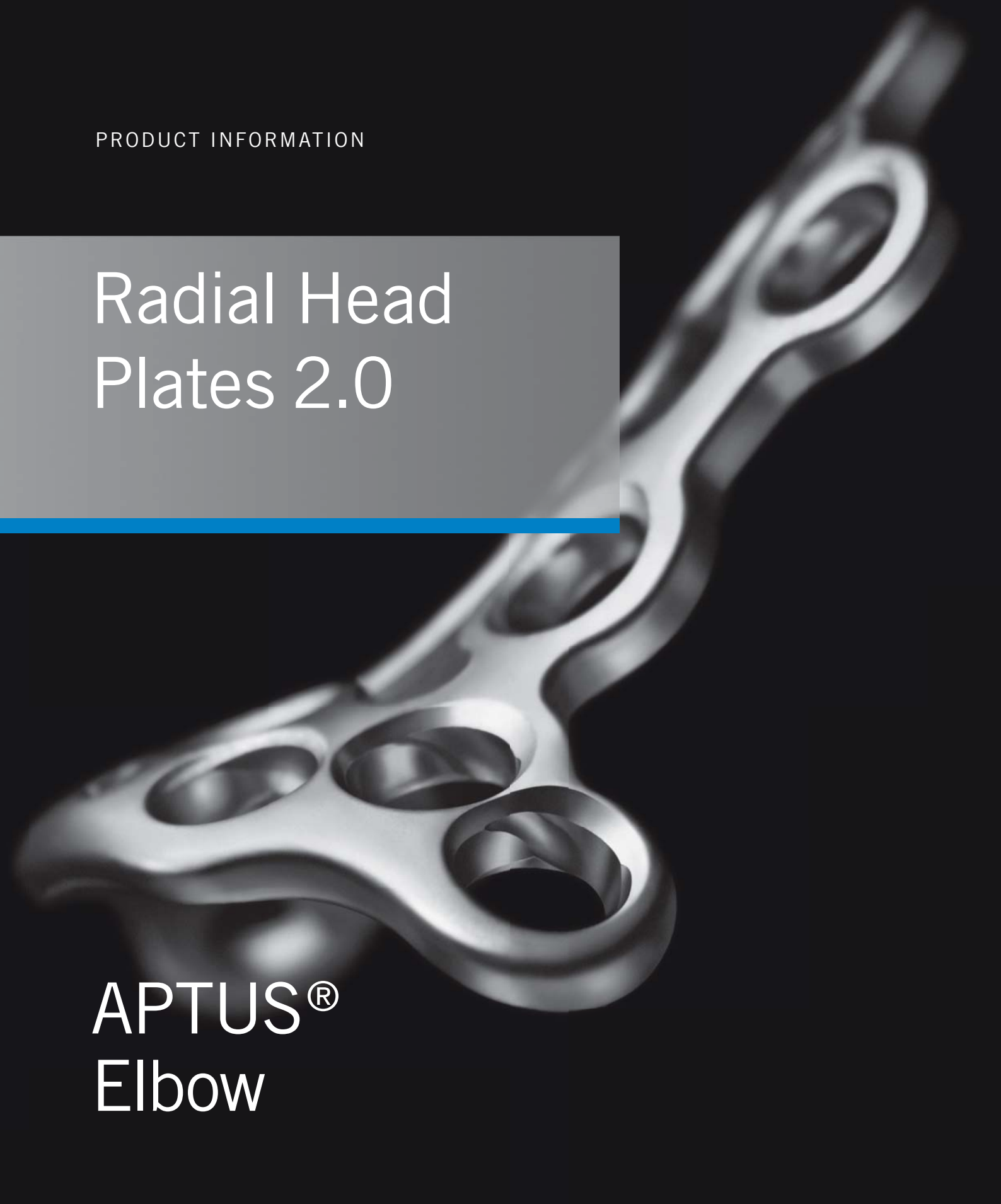


# Radial Head Plates 2.0

Safe and stable internal fixation of complex radial head fractures with anatomic angular stable plates

## Plate Features

- Minimal plate thickness
- Chamfered plate contour to minimize soft tissue irritation
- Anatomic plate design
- TriLock – multidirectional ( $\pm 15^\circ$ ) and angular stable screw locking is possible in each plate hole
- All plate holes can be utilized with either locking or non-locking screws
- Consistent screw size ( $\varnothing$  2.0 mm)
- Titanium Grade 4 for optimal stability
- Holes for temporary 1.2 mm K-wire fixation

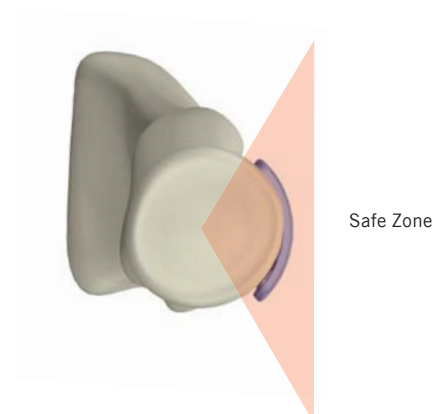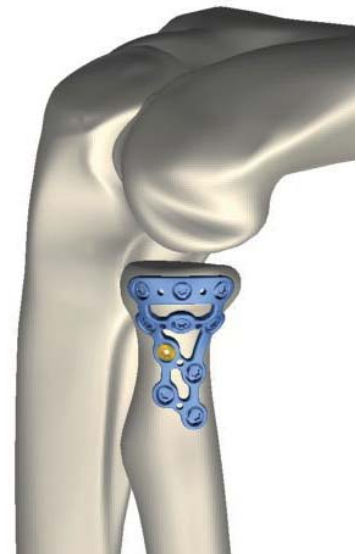

## Clinical Benefits

### Anatomic design to minimize soft tissue irritation and ensure stable fixation

- Anatomically optimized plate designs
- Fixation of multiple fragments
- Low plate profile minimizes soft tissue irritation, particularly in the case of the annular ligament



**Radial Head Rim Plate**  
For the stabilization of complex comminuted fractures

### Multiple locking options

- Multidirectional and angular stable TriLock technology provides a variable cone of  $\pm 15^\circ$  for optimal fixation of all fragments
- Angular stability of screw fixation minimizes risk of post-operative loss of reduction
- Screw configuration can be adapted to the fracture pattern due to numerous plate holes



**Radial Head Buttress Plate**  
For an optimal sparing of the annular ligament


## Clinical Cases

### Case 1



Radial Head Rim Plate for a Mason Typ III / IV head fracture

### Case 2

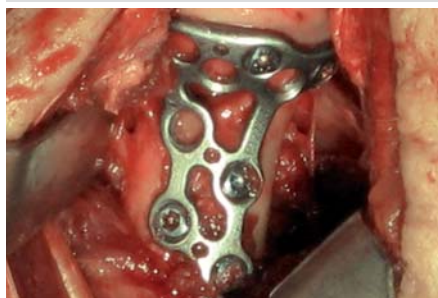


Radial Head Buttress Plate for a Mason Type III neck fracture

# Technology, Biomechanics, Screw Features


## Multidirectional and angular stable TriLock locking technology

### Technology

- Multidirectional ( $\pm 15^\circ$ ) and angular stable TriLock locking system
  - Spherical three-point wedge-locking
  - Friction locking through radial bracing of the screw head in the plate – without additional tensioning components
- TriLock screws can be re-locked in the same plate hole under individual angles up to three times
- Minimal screw head protrusion thanks to internal locking contour
- No cold welding between plate and screws
- Intra-operative fine tuning capabilities

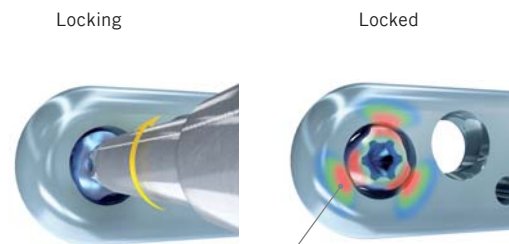

Secure locking of the TriLock screw

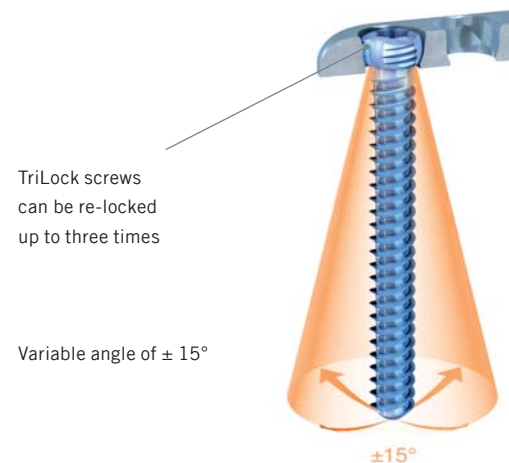


TriLock screws can be re-locked up to three times

Variable angle of  $\pm 15^\circ$


$\pm 15^\circ$

Minimal screw head protrusion

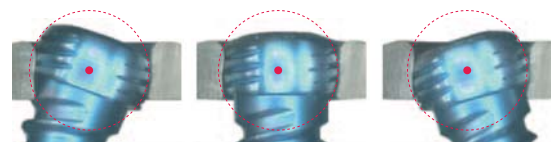



## Biomechanics

- Internal fixator principle
  - Forces around the proximal radius bypass the generally unstable fracture site
  - Low contact for ideal blood supply
  - Functionally dynamic construct to avoid cut-out of screws

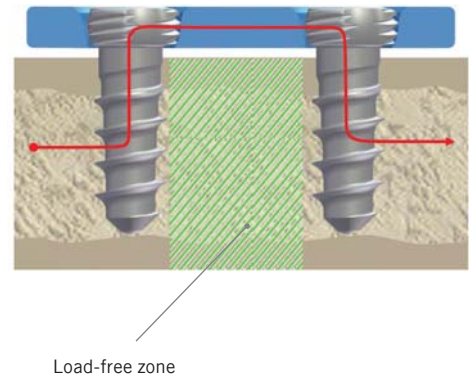

## Screw Features

- HexaDrive screw head design
  - Secure connection between screw and screwdriver
  - Increased torque transmission
  - Optimal self-retaining mechanism
- Maximum soft tissue protection due to chamfered shape of the screw head without sharp edges
- Atraumatic tip prevents soft tissue irritation when inserting screws bicortically
- Tapered core diameter for increased torsional and tensile strength
- Precision cut thread profile for improved sharpness and self-tapping properties
- Double threaded for faster insertion of TriLock screws
- TiAl6V4 for improved strength

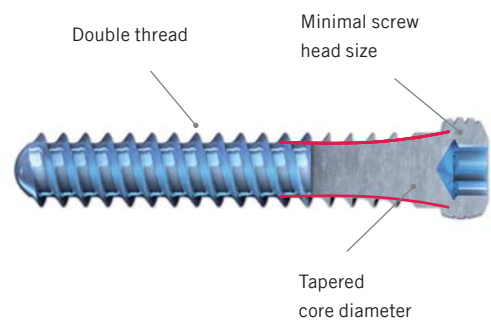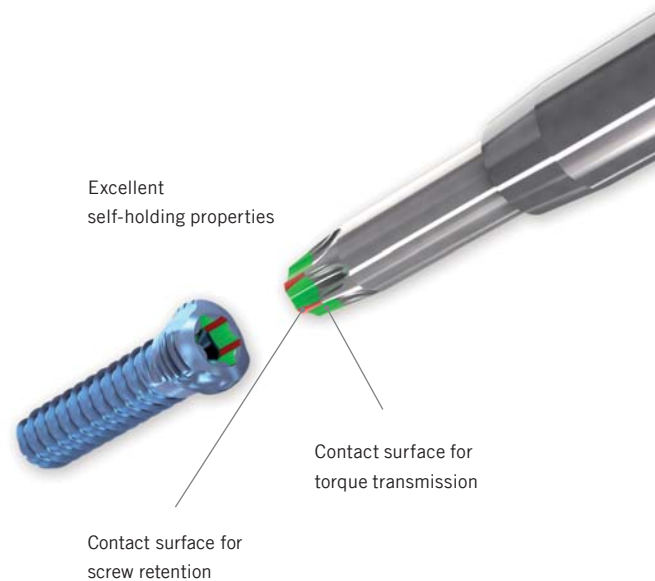

# Storage in Perfection

- Modular concept
- Economic and compact system
- Easy to handle
- Streamlined organization of implants and instruments



APTUS Radial Head Plates 2.0  
Implant Container (incl. Stickers, empty)


Art. No. A-0450.01



# Ordering Information

## 2.0 TriLock Radial Head Plates

Material: Titanium (ASTM F67)  
Plate thickness: 1.4 mm

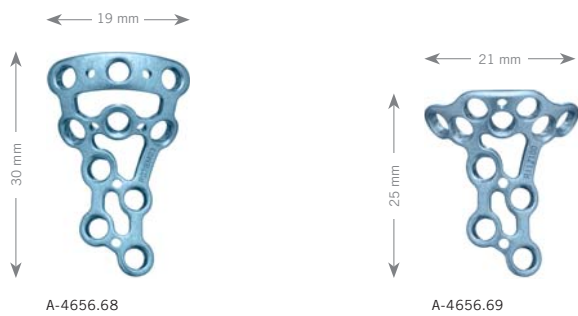

Art. No.	Description	Holes	Piece/Pkg
A-4656.68	rim plate	10	1
A-4656.69	buttress plate	11	1

## K-Wires, Stainless Steel



Art. No.	System Size	Description	Length	Piece/Pkg
A-5040.21	1.2	trocar	150 mm	10
A-5042.21	1.2	lancet	150 mm	10


## K-Wire Dispenser



Art. No.	System Size	Length	Piece/Pkg
A-6010.12	1.2	185 mm	1

## 2.0 TriLock Screws, HexaDrive 6

Material: Titanium (ASTM F136)

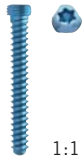

1:1

Length	Art. No.	Piece/Pkg	Art. No.	Piece/Pkg
6 mm	A-5450.06/1	1	A-5450.06	5
8 mm	A-5450.08/1	1	A-5450.08	5
10 mm	A-5450.10/1	1	A-5450.10	5
12 mm	A-5450.12/1	1	A-5450.12	5
14 mm	A-5450.14/1	1	A-5450.14	5
16 mm	A-5450.16/1	1	A-5450.16	5
18 mm	A-5450.18/1	1	A-5450.18	5
20 mm	A-5450.20/1	1	A-5450.20	5
22 mm	A-5450.22/1	1	A-5450.22	5
24 mm	A-5450.24/1	1	A-5450.24	5
26 mm	A-5450.26/1	1	A-5450.26	5
28 mm	A-5450.28/1	1	A-5450.28	5
30 mm	A-5450.30/1	1	A-5450.30	5

## 2.0 Cortical Screws, HexaDrive 6

Material: Titanium (ASTM F136)



1:1

Length	Art. No.	Piece/Pkg	Art. No.	Piece/Pkg
6 mm	A-5400.06/1	1	A-5400.06	5
8 mm	A-5400.08/1	1	A-5400.08	5
10 mm	A-5400.10/1	1	A-5400.10	5
12 mm	A-5400.12/1	1	A-5400.12	5
14 mm	A-5400.14/1	1	A-5400.14	5
16 mm	A-5400.16/1	1	A-5400.16	5
18 mm	A-5400.18/1	1	A-5400.18	5
20 mm	A-5400.20/1	1	A-5400.20	5
22 mm	A-5400.22/1	1	A-5400.22	5
24 mm	A-5400.24/1	1	A-5400.24	5
26 mm	A-5400.26/1	1	A-5400.26	5
28 mm	A-5400.28/1	1	A-5400.28	5
30 mm	A-5400.30/1	1	A-5400.30	5

## 2.3 Cortical Screws, HexaDrive 6

Material: Titanium (ASTM F136)



1:1

Length	Art. No.	Piece/Pkg	Art. No.	Piece/Pkg
8 mm	A-5500.08/1	1	A-5500.08	5
10 mm	A-5500.10/1	1	A-5500.10	5
12 mm	A-5500.12/1	1	A-5500.12	5
14 mm	A-5500.14/1	1	A-5500.14	5
16 mm	A-5500.16/1	1	A-5500.16	5

## Twist Drills Ø 1.6 mm



A-3414



A-3424



A-3434

Art. No.	System Size	Stop	Length	Drill Shaft End	Piece/Pkg
A-3414	2.0	30 mm	92 mm	Dental	1
A-3424	2.0	30 mm	92 mm	Stryker J-Latch	1
A-3434	2.0	30 mm	86 mm	AO Quick Coupling	1

## Twist Drills Ø 2.1 mm (for gliding holes)



A-3411



A-3421



A-3431

Art. No.	System Size	Stop	Length	Drill Shaft End	Piece/Pkg
A-3411	2.0	10 mm	72 mm	Dental	1
A-3421	2.0	10 mm	72 mm	Stryker J-Latch	1
A-3431	2.0	10 mm	66 mm	AO Quick Coupling	1

## Countersink for Cortical Screws



Art. No.	System Size	Ø	Length	Shaft End	Piece/Pkg
A-3610	2.0/2.3	3.5	37 mm	Dental	1

## Handle with Quick Connector



Art. No.	Length	Shaft End	Piece/Pkg
A-2071	107 mm	Dental	1

## Drill Guide



Art. No.	System Size	Description	Length	Piece/Pkg
A-2020	1.2 - 2.3	centric, excentric	149 mm	1

## Depth Gauge



Art. No.	System Size	Length	Piece/Pkg
A-2032	2.0 / 2.3	151 mm	1

## Screwdriver, self-holding



Art. No.	System Size	Interface	Length	Piece/Pkg
A-2610	2.0 / 2.3	HD6	153 mm	1

## Plate Bending Pliers



Art. No.	System Size	Description	Length	Piece/Pkg
A-2040	1.2 - 2.3	with Vario pin	119 mm	1

### Instrument Trays



A-6020



A-6025

Art. No.	Description	Piece/Pkg
A-6020	with silicon inlay	1
A-6025	for A-6020	1

### Lid for Instrument Tray



Art. No.	Description	Piece/Pkg
M-6198	for A-6020	1

ELBOW-01000001\_v3 / © 06.2011, Medartis AG, Switzerland. All technical data subject to alteration.

#### **HEADQUARTERS**

Medartis AG | Hochbergerstrasse 60E | 4057 Basel/Switzerland  
P +41 61 633 34 34 | F +41 61 633 34 00 | [www.medartis.com](http://www.medartis.com)

#### **SUBSIDIARIES**

Australia | Austria | France | Germany | Mexico | New Zealand | Poland | UK | USA

For detailed information regarding our subsidiaries and distributors, please visit [www.medartis.com](http://www.medartis.com).