


medartis

PRECISION IN FIXATION

PRODUCT INFORMATION

# All-in-One Staple

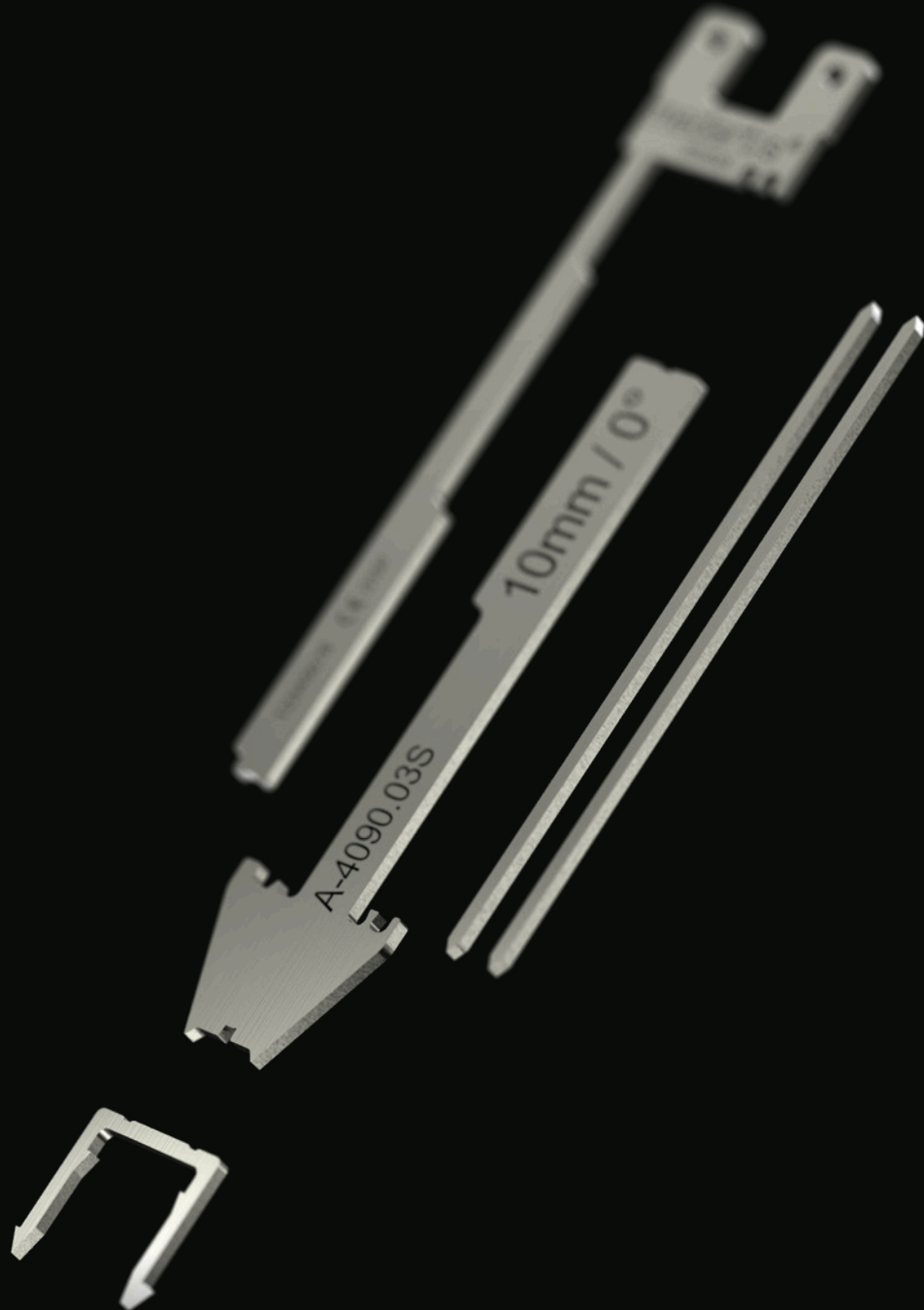

**APTUS** Foot


# Ready. Steady. Staple.

The pocket size sterile staple kit with implant and all implant-specific instruments in a lightweight packaging of 19 g and 66 x 99 x 11 mm.

A further step towards comprehensive hallux treatment and bone fixation in the foot.


## The Swiss Pocket Staple

- Complete
- Functional
- Sterile and compact on the shelf

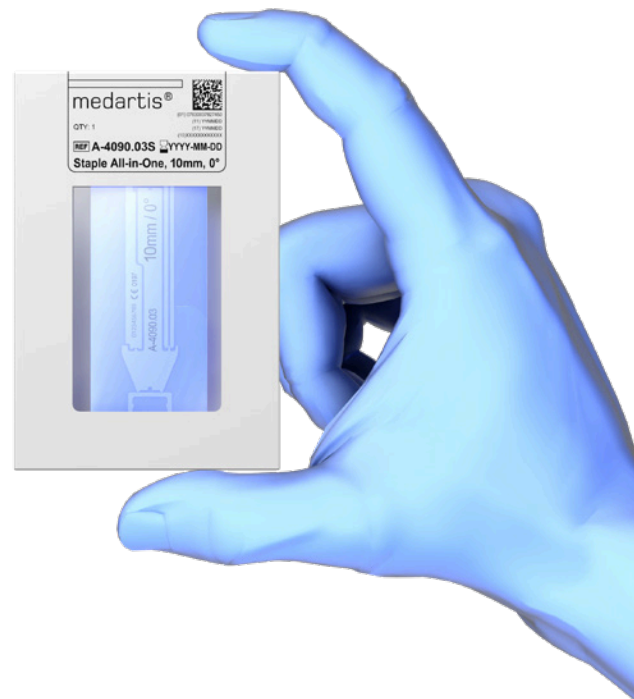

## Implant and Instruments

- 1-piece precision laser cut from stainless steel sheet material
- Implant and instruments can be easily separated
- Different staple designs to meet anatomical needs and surgeon's preference
- Swiss engineering and manufacturing

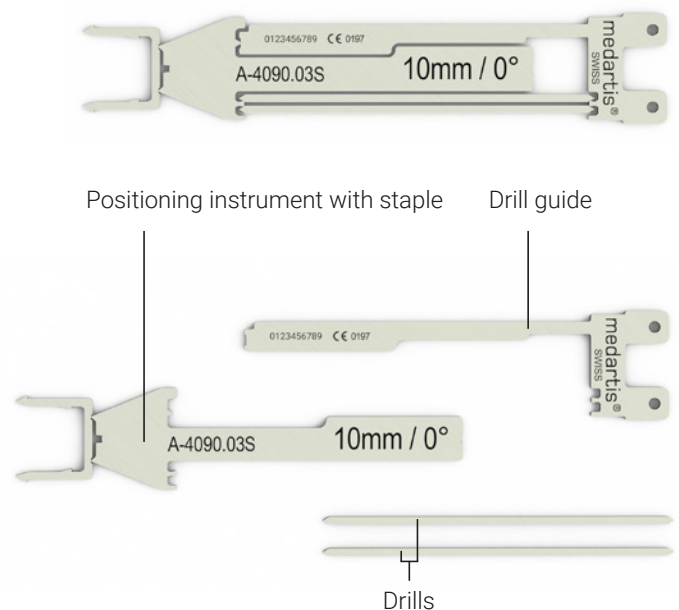

## Specifications

- Available in two widths (8 mm and 10 mm) and two angles (0° and 26°)
- Leg length 10 mm
- Barbed legs can prevent migration of the staple
- Staple design and drill guide allow for compression
- Two drills and the corresponding drill guide are included
- Positioning instrument for precise insertion



Staple All-in-One  
8 mm, 0°



Staple All-in-One  
8 mm, 26°



Staple All-in-One  
10 mm, 0°



Staple All-in-One  
10 mm, 26°

## Values

- Saving reprocessing capacity and resources
- Lean material and supply management – 1 item, 1 price
- Reducing OR idle time

R FOOT-07000001\_v0 / © 2021-03, Medartis AG, Switzerland. All technical data subject to alteration.

#### **MANUFACTURER & HEADQUARTERS**

Medartis AG | Hochbergerstrasse 60E | 4057 Basel / Switzerland  
P +41 61 633 34 34 | F +41 61 633 34 00 | [www.medartis.com](http://www.medartis.com)

#### **SUBSIDIARIES**

Australia | Austria | Brazil | China | France | Germany | Japan | Mexico | New Zealand | Poland | Spain | UK | USA

For detailed information regarding our subsidiaries and distributors, please visit [www.medartis.com](http://www.medartis.com)



Disclaimer: This information is intended to demonstrate the Medartis portfolio of medical devices. A surgeon must always rely on her or his own professional clinical judgement when deciding whether to use a particular product when treating a particular patient. Medartis is not giving any medical advice. The devices may not be available in all countries due to registration and / or medical practices. For further questions, please contact your Medartis representative ([www.medartis.com](http://www.medartis.com)). This information contains CE-marked products.  
For US only: Federal law restricts this device to sale by or on the order of a physician.